

4040XL Heavy duty examination table

An innovative general examination table for demanding use



Product pictured with optional accessories.



Main features:

1. The patient surface width (up to 100cm) together with a thick upholstery ensure a pleasant experience
2. A very high safe working load (SWL) of 300kg and new features improve patient safety
3. Height adjustment range of 42 – 97 cm helps patients and makes working more ergonomic

The new 4040XL is an excellent choice for all healthcare facilities that need to treat also overweight and obese patients. Its patented construction and extremely low minimum height ensure that even patients with wheelchair or obese patients are able to get up onto the table with ease. High maximum height improves ergonomics of the care providers. The patient surface width together with a thick upholstery (foam 6 cm) provide very pleasant experience for the patient. The default patient surface width is 90 cm (80 cm and 100 cm as options).

In all Lojer products patient safety is of utmost importance. The 4040XL has a new safety feature: the table top rises automatically if there's something caught between the mechanism. The very high safe working load (SWL) of 300 kg makes it safe to treat also bariatric patients in all positions and heights.


The 4040XL examination table has a new digital control unit, which notifies user of maintenance need. It also allows locking up the table with the hand control. The table features a paper roll holder and cutter as standard, as well as an adapter for fitting an IV-stand (6 possible fitting positions). A two-piece hands-free height adjustment bar is available as an option, enabling the table height and back section to be adjusted ergonomically and hygienically from both sides of the table by foot. Also electric Trendelenburg and anti-Trendelenburg are available as options.

LOJER® For easy care

Technical data and features	
Total width	95 cm with 90 cm upholstery 80 cm with 80 cm upholstery 100 cm with 100 cm upholstery
Patient surface width	90 cm (80 / 100 cm optional)
Total length	200 cm
Patient surface length	200 cm
Total weight	120 kg
Height adjustment	42 - 97 cm
Back section length	70 cm
Back section adjustment	73°
Electric Trendelenburg and anti-Trendelenburg adjustment (optional)	max. 20°
Castors	Ø125 mm (125/150 mm twin castors optional)
Safe working load (SWL)	300 kg
Upholstery	Upholstery thickness 75 mm (foam 60 mm)
Other features	Hand control paper roll holder and cutter for 80 cm paper adapter for fitting an IV-stand (6 possible fitting positions)
CE -marking	Yes
Warranty	2/10 years
Country of Manufacture	Finland

Factory installed options:

- Two-piece hands-free adjustment bar for height and back section
- Electric Trendelenburg and anti-Trendelenburg (battery included)
- Battery power, only for emergency use
- Central locking twin castors Ø125mm / Ø150mm
- Patient surface width 80cm and 100cm

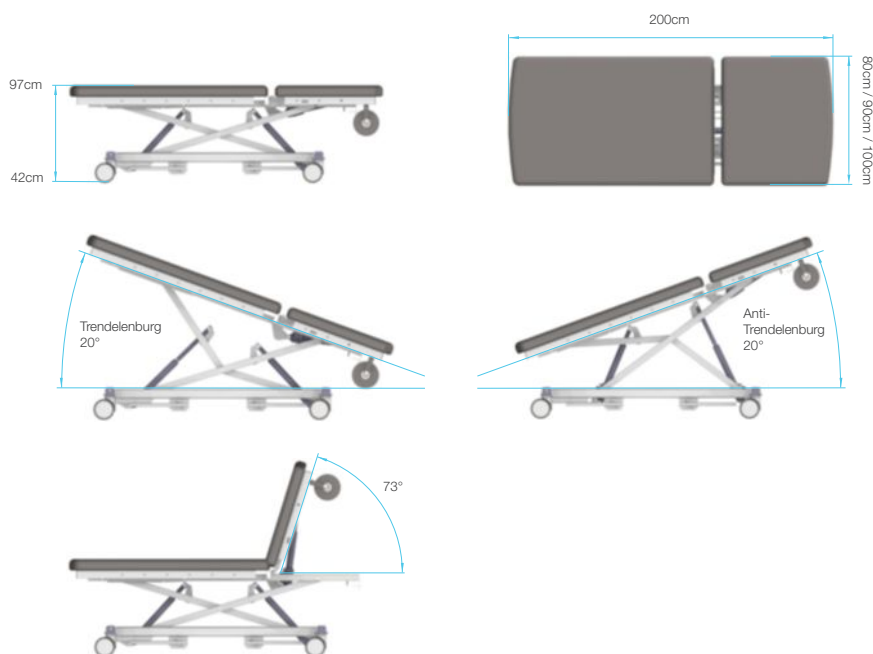

Other accessories:

- Side rails
- Foot control (FS32)
- Accessory rail
- Examination light
- IV-pole
- Neck cushion
- Push handles
- Anaesthesia frame
- Arm support
- Adapter for accessories
- Clamps
- All accessories fitted to accessory rail



LOJER® For easy care